


LIS

The Phoenix M50 can be connected directly to a Laboratory Information System (LIS). Note that LIS is a general acronym; it is possible that a system might also be referred to as a HIS (Hospital Information System).

In a normal LIS setup, there is a LIS server somewhere in the facility that is connected to a terminal server. BD Instruments connect to the terminal server via serial line. The facility may use RJ45 and CAT 5 cabling, but the communications protocol should be serial for the BD LIS interface to work. In the case of a RJ45 connector, a RJ45 to DB9 Serial adapter will be required. If the facility uses a standard RS232 DB9 connector, then it can be plugged in directly.

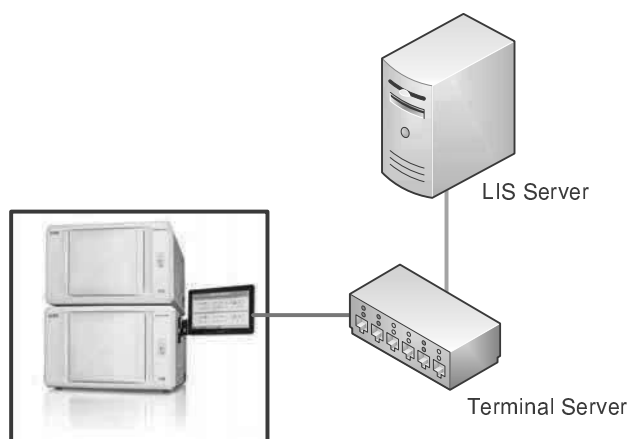


Connections

The LIS connects to the DB9 Serial Port on the bottom of the All In One Computer.




This is a standard RS-232 DB-9 connector with standard pinouts. The Phoenix M50 uses Pin 2 (Receive), Pin 3 (Transmit), and Pin 5 (Ground). If the facility uses CAT 5 Cabling for this connection, an adapter may be required.

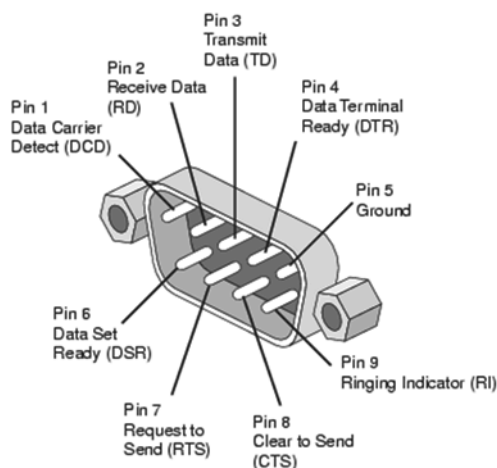

User Software Setup

From the **Configuration -> Communications Configuration** tabs check LIS Enabled and **Save**.

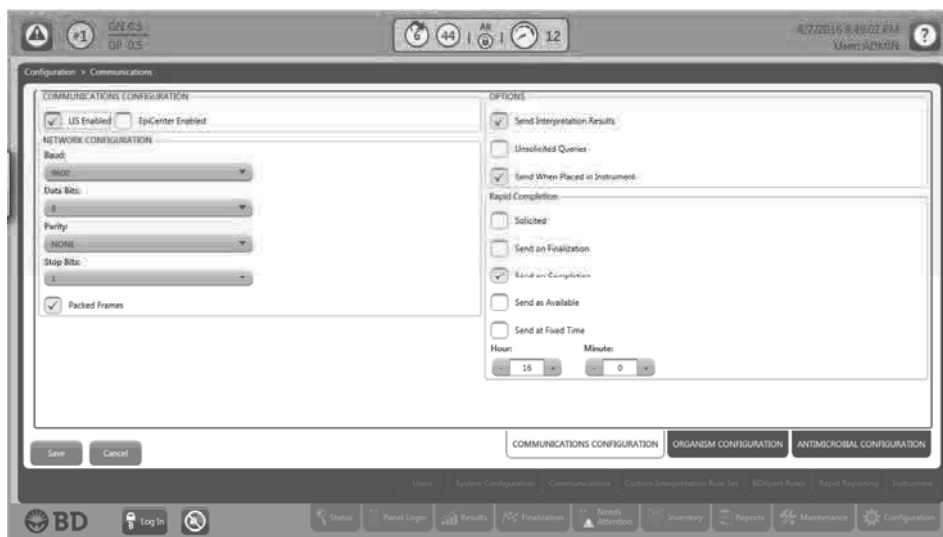

The low level LIS options are now available. These are standard and include Baud Rate, Data Bits, Stop Bits, Parity, and Packed/Unpacked Frames.

Refer to the LIS Vender specification for more detailed description of the BD LIS Interface.