


Laboratory Information System (LIS) Connection (for Getein1100)

1. Serial Port Driver Installation

- a) Prepare RS232 data line and driver software, as the following Figure.



- b) Put the above disc into CD-ROM to drive correctly and check the detail information.
- c) Double-click “ autorun” to install RS232 driver software automatically.

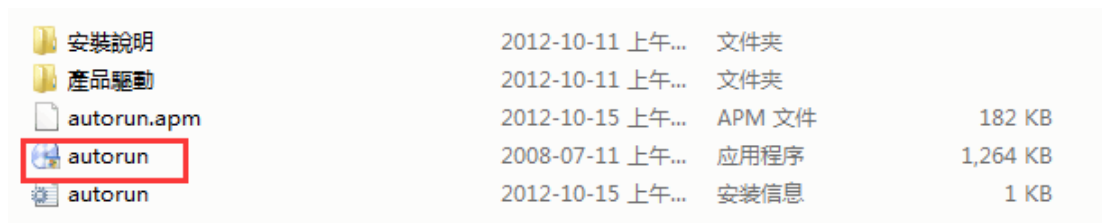

2. Laboratory Information System Connection

- a) Turn on Getein1100, and connect the analyzer and computer through RS232 data line.




- b) Keep Communication Status “ON” and select “Serial port” in communication mode.

No.9 Bofu Road, Luhe District, Nanjing (211505) China.

Tel: 86-25-68568707

Fax: 86-25-68568500

Website: www.bio-GP.com.cn

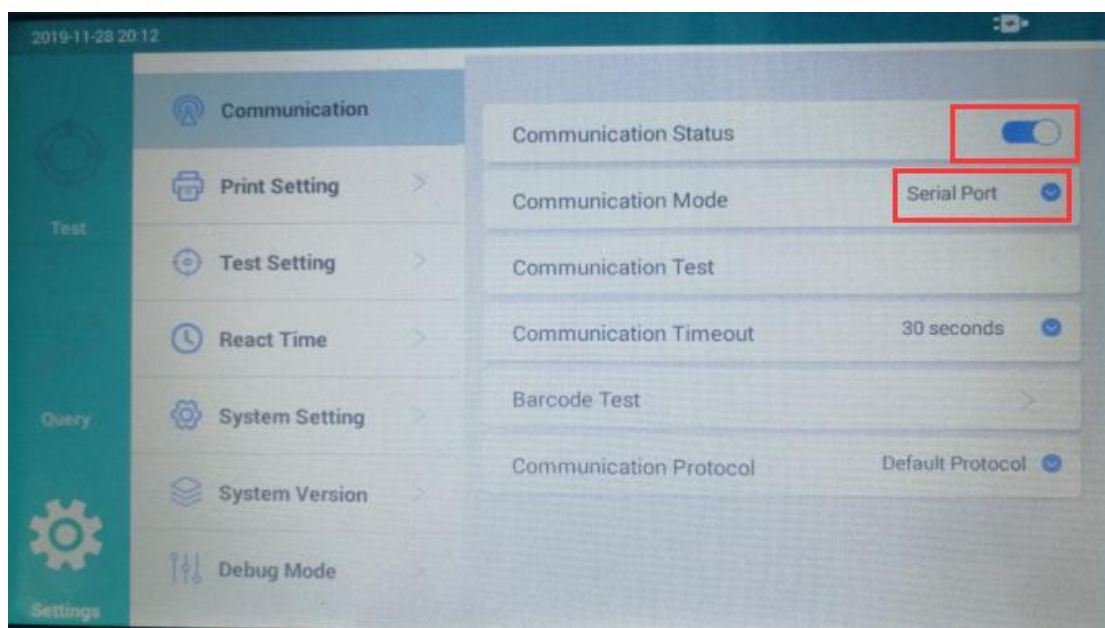

c) Select the communication protocol as needed (Getein can provide the communication protocols, hospital LIS providers help the LIS connection).

