

YHLO SHENZHEN YHLO BIOTECH CO., LTD.

VISION Pro LIS Configuration manual

V1.0

Contents

1. Overview.....	1
1.1. Introduction.....	1
1.2. Manual Scope.....	1
2. LIS Connection Configuration.....	2
2.1. LIS Server IP Address & Port Setting.....	2
2.2. Device IP Setting.....	2
2.3. Other settings.....	3
2.4. LIS Connection Confirm.....	3
2.5. Notes.....	3

1. Overview

1.1. Introduction

Through the correct configuration of LIS, the VISION ESR Analyzer can automatically send sample results and quality control results to the LIS server in real time or manually, individually or in batches.

1.2. Manual Scope

This manual is only applicable to VISION ESR Analyzer LIS configuration.

2. LIS Connection Configuration

2.1. LIS Server IP Address & Port Setting

- 1) Log in Admin account and enter 'Settings - SysSet - LIS' interface as showed below.

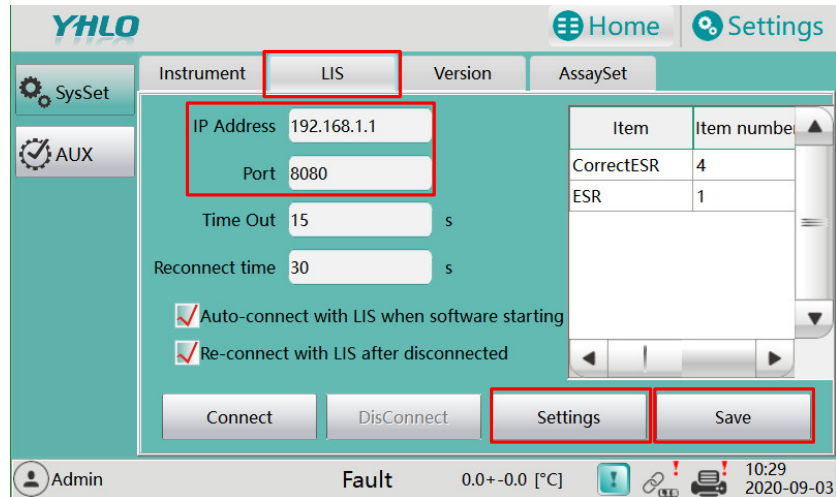


Fig. 1 LIS Server IP Address & Port Setting

- 2) Input IP Address and Port of LIS server, then click Save button.

2.2. Device IP Setting

- 1) Click Settings button in Fig. 1, enter device IP setting interface.
- 2) Input device IP address, Sub Net Mask and Gateway, then click Save. (You may ask the network administrator for IP distributing.)

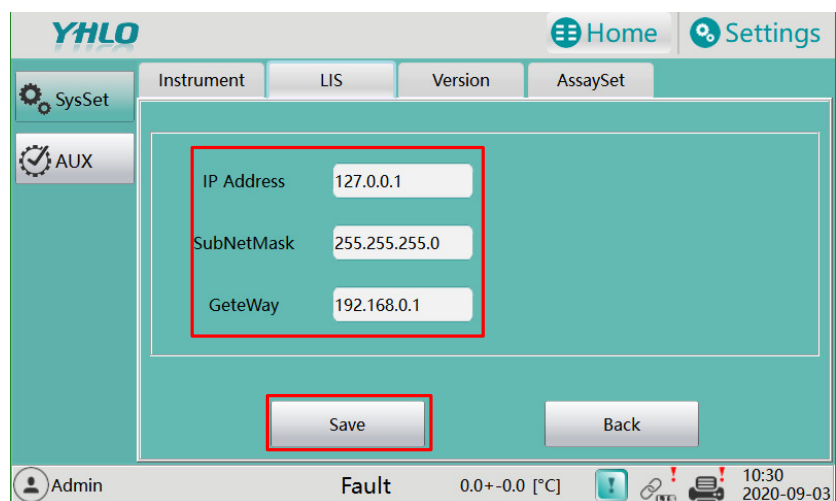


Fig. 1 Device IP Setting

2.3. Other settings

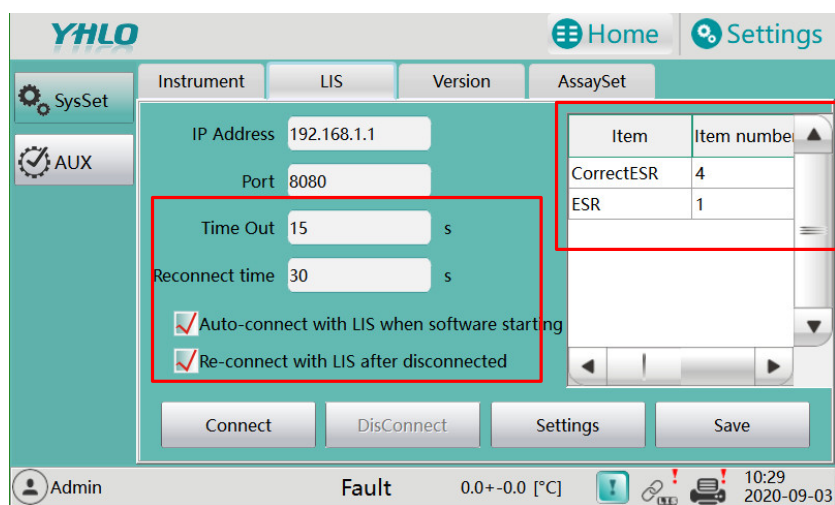


Fig.3 Other Settings

- 1) Time Out: Set the timeout limit of the connection attempt between the instrument and LIS, the input range is 10~60s.
- 2) Reconnect time: Set the interval time for the system to automatically reconnect to the LIS after the LIS connection is disconnected. The input range is 10~60s.
- 3) Auto-connect with LIS when software starting: After checking this option, the system will establish a connection with the LIS server when the system starts.
- 4) Re-connect with LIS after disconnected: After checking this option, the system will check the LIS connection status according to 'Reconnect time'. If LIS is disconnected, it will automatically reconnect to the LIS.
- 5) Item number: Item and Item Number are fixed, CorrectESR represents result with temperature correction, ESR represents result without temperature correction.

2.4. LIS Connection Confirm

After configuration, click the Connect button. If the connection is successful, the exclamation mark of the LIS icon will disappear, as shown in the figure below.



2.5. Notes

- 1) Before setting up LIS, click the Disconnect button to disconnect the LIS first.
- 2) After setting, click the Save button.