

XN-L Interfacing Guide

HCLAB WAM

Author: Lawrence Tay

Last Edited: 23 January 2015

Version: 1.0

Document Properties

Document information

Document type	Manual
Version	1.0
Status	Draft
Last edited	23-Jan-15
Author	Lawrence Tay
File Name	XN-L Interfacing Guide

Document history

Version	Date	Author	Comments
1.0	01-Jan-15	Lawrence Tay	

Document distribution

Sysmex Asia Pacific Pte Ltd

Contents

Pre-Requisites	4
Setup Steps	5
Configure Settings in Setting	5
Host Computer settings	5
QC Auto Output Setting	9
Operational Setting	11

Pre-Requisites

Before enabling interfacing settings , please ensure the following conditions are met:

1. You should obtain Service credentials for login IPU software.
2. Installation of IPU software is completed.

Setup Steps

The basic setup steps are as follows:

1. Configuration changes under “Setting” button.
 - a. Host Computer settings
 - b. QC auto output settings
 - c. Operational Setting
 - i. Auto Validate
 - ii. Analysis Ordering
 - iii. Auto Output
2. Configuration Settings in the IPU Settings
3. Configuration Settings in the IPU Analyser Settings


Configure Settings in Setting

Host Computer settings


1. Log in the service account in the IPU software.



2. Go to main Menu screen.

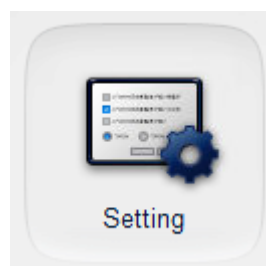


3. Click on the Setting button.

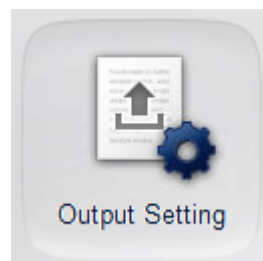


4. Click Output Setting button.


5. Click Host Computer button.

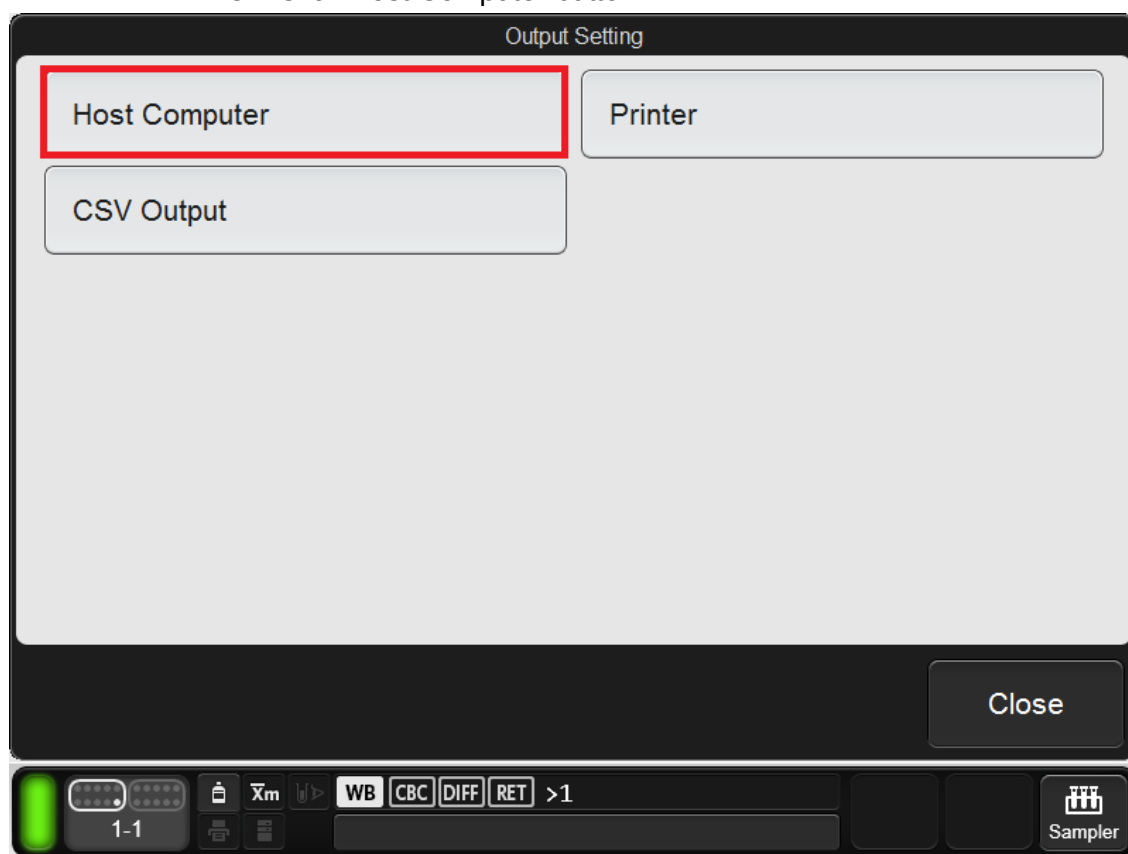

6. Connect to Host Computer using TCP/IP communication,
 - a. Tick - Connect to Host Computer
 - b. Select - TCP/IP Connection
 - c. Click - Modify Settings button

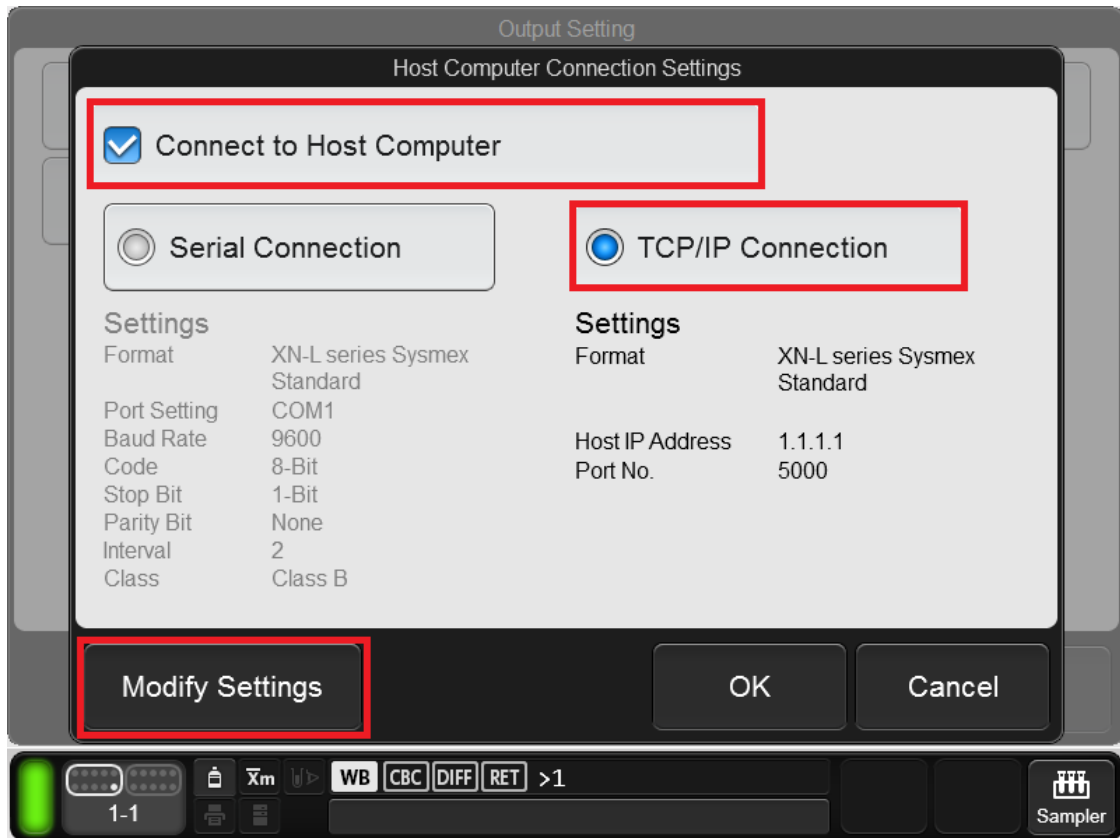

7. TCP/IP Connection settings,
 - a. Format : XN-L Series ASTM1381-02/ASTM1894-97
 - b. Host IP Address : e.g. 192.168.28.200
(IP Address of HCLAB WAM/LIS)
 - c. Port No. : e.g. 5000
(Port No. used for communication between XN-L and HCLAB WAM)
 - d. Click on OK button to save the changes.

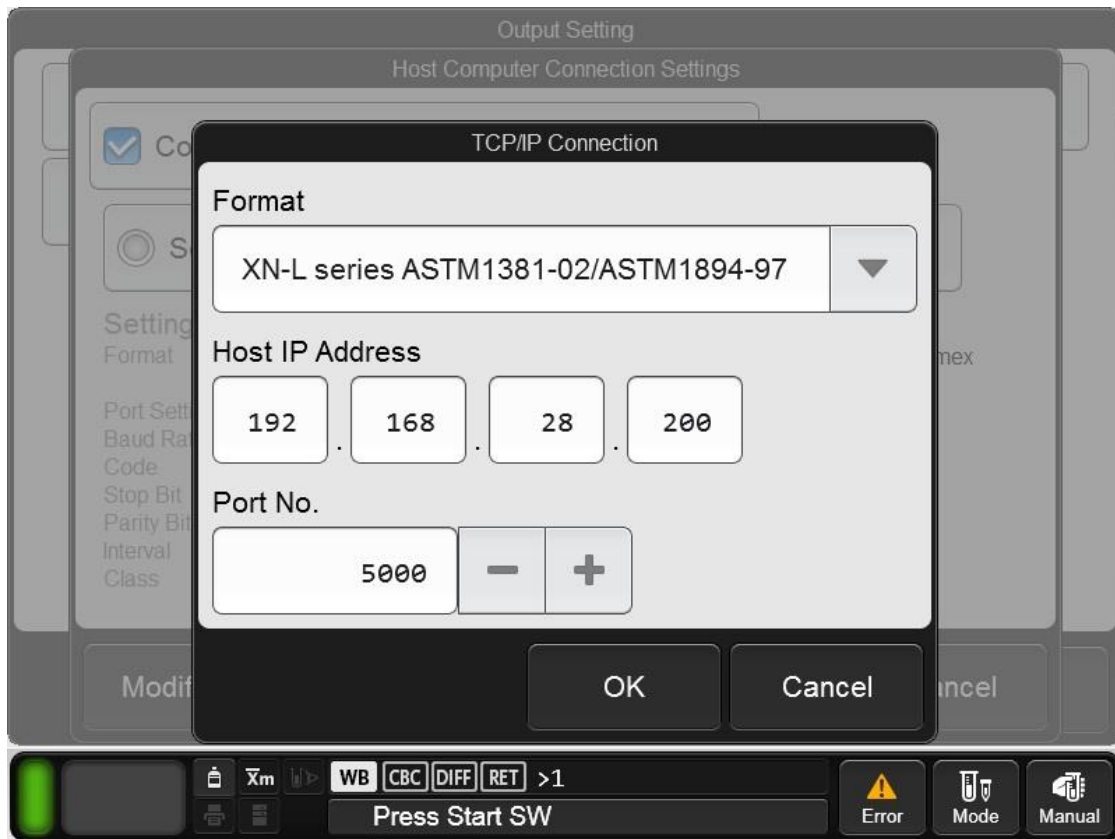


Note: You should do a ping test to check that network communication is working.

QC Auto Output Setting



1. Go back to main Menu screen.

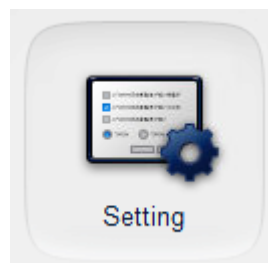


2. Click on the **"Setting"** button.




3. Click on **"QC Setting"** button.

4. Click on **"QC Chart Data Auto Output Setting"** button

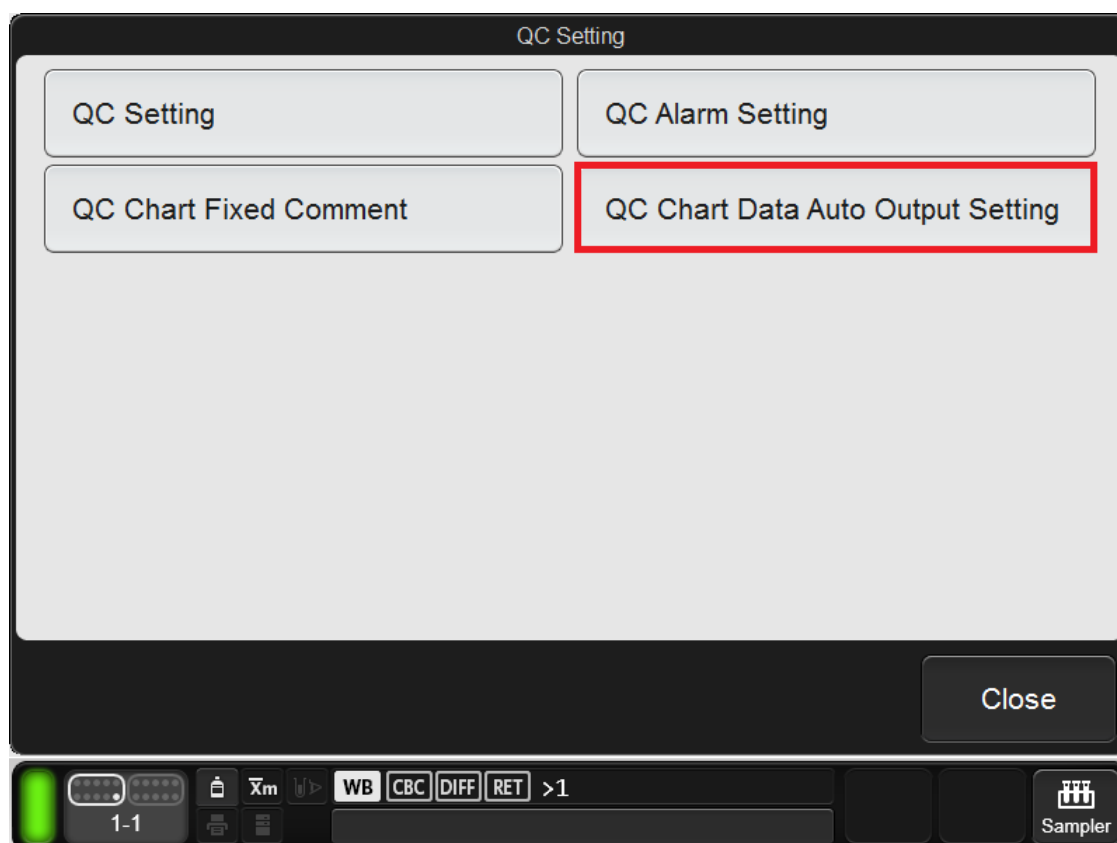

5. Ensure the QC settings as shown below,
 - a. Tick – X-barM Files
 - b. Tick – Host Computer (HC)
 - c. Click on OK button to save the changes.

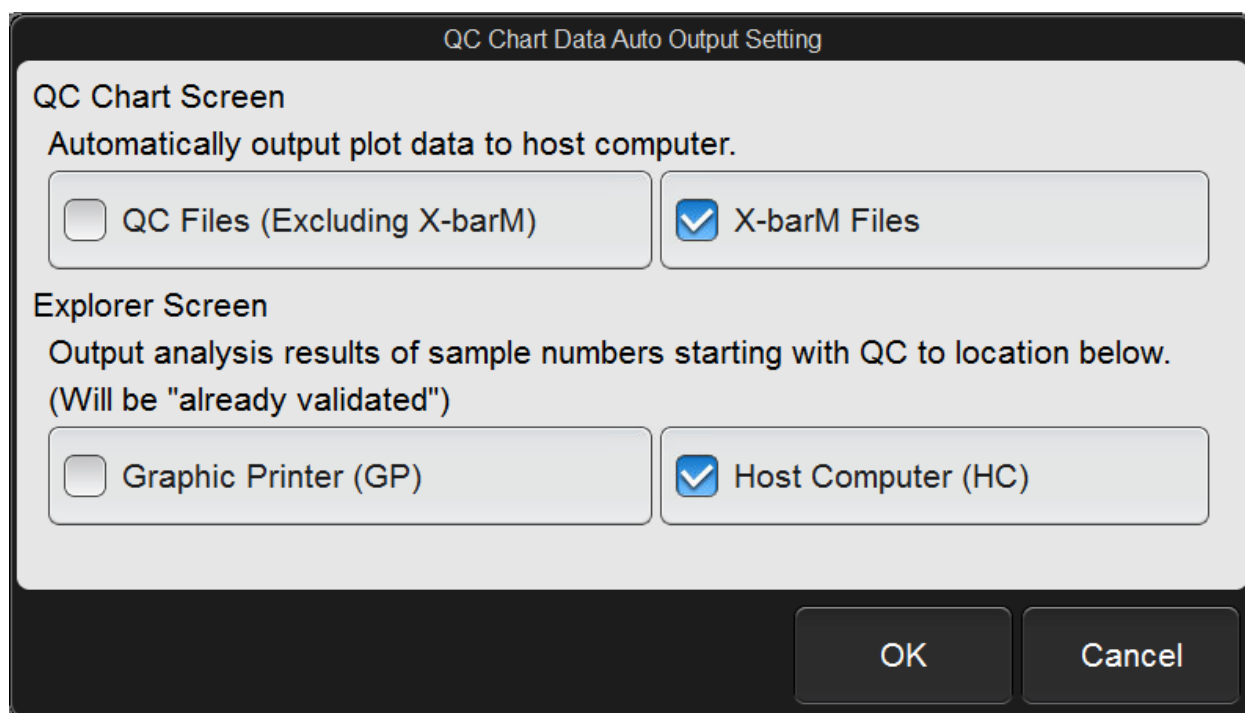

The image shows a dialog box titled "QC Chart Data Auto Output Setting". It contains two sections: "QC Chart Screen" and "Explorer Screen".

QC Chart Screen
Automatically output plot data to host computer.

There are two checkboxes in this section:
☐ QC Files (Excluding X-barM)
☒ X-barM Files

Explorer Screen
Output analysis results of sample numbers starting with QC to location below.
(Will be "already validated")

There are two checkboxes in this section:
☐ Graphic Printer (GP)
☒ Host Computer (HC)


At the bottom right, there are two buttons: "OK" and "Cancel".

Operational Setting


1. Go back to main Menu screen.




2. Click on the Setting button.

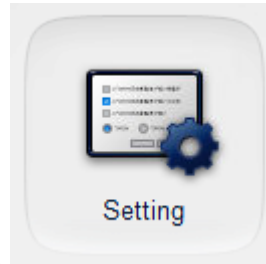


3. Click on Operational Settings button.

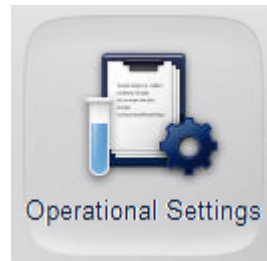


4. Click – Auto Validate button

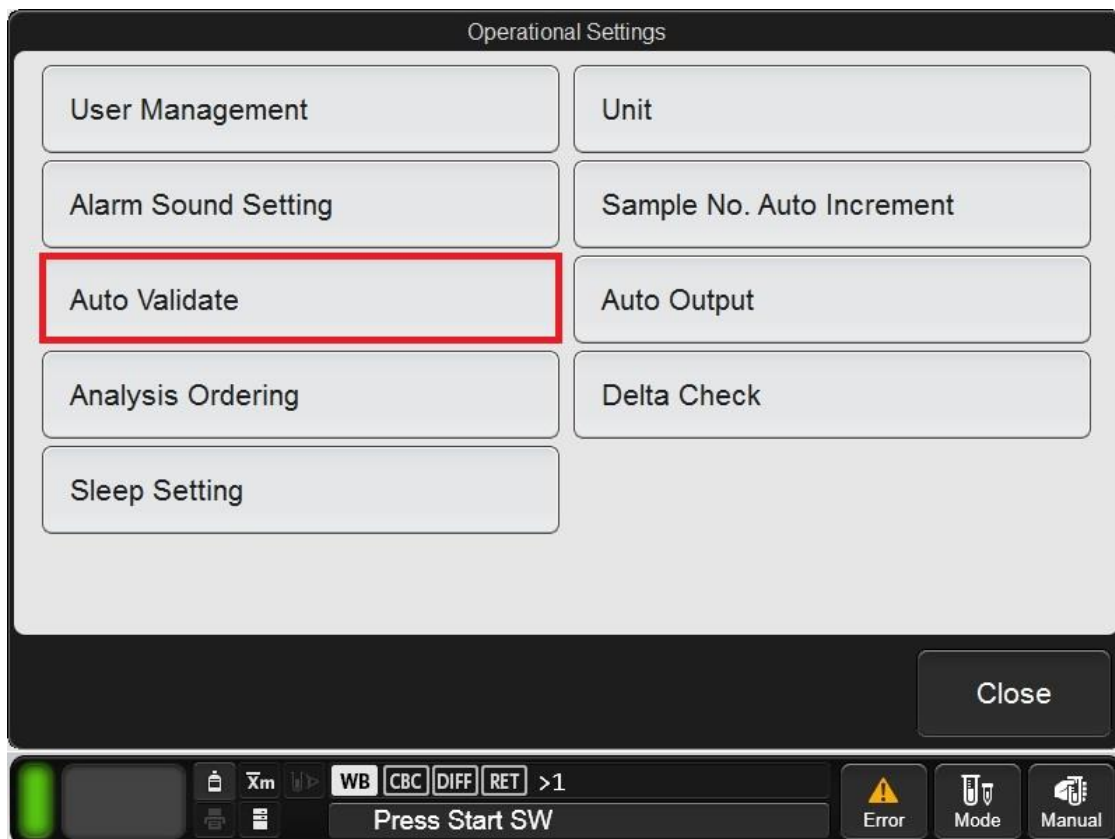


5. Set the settings as shown below
 - a. Tick – Auto Validate
 - b. Select – Use simple settings
 - c. Click - Modify Settings button

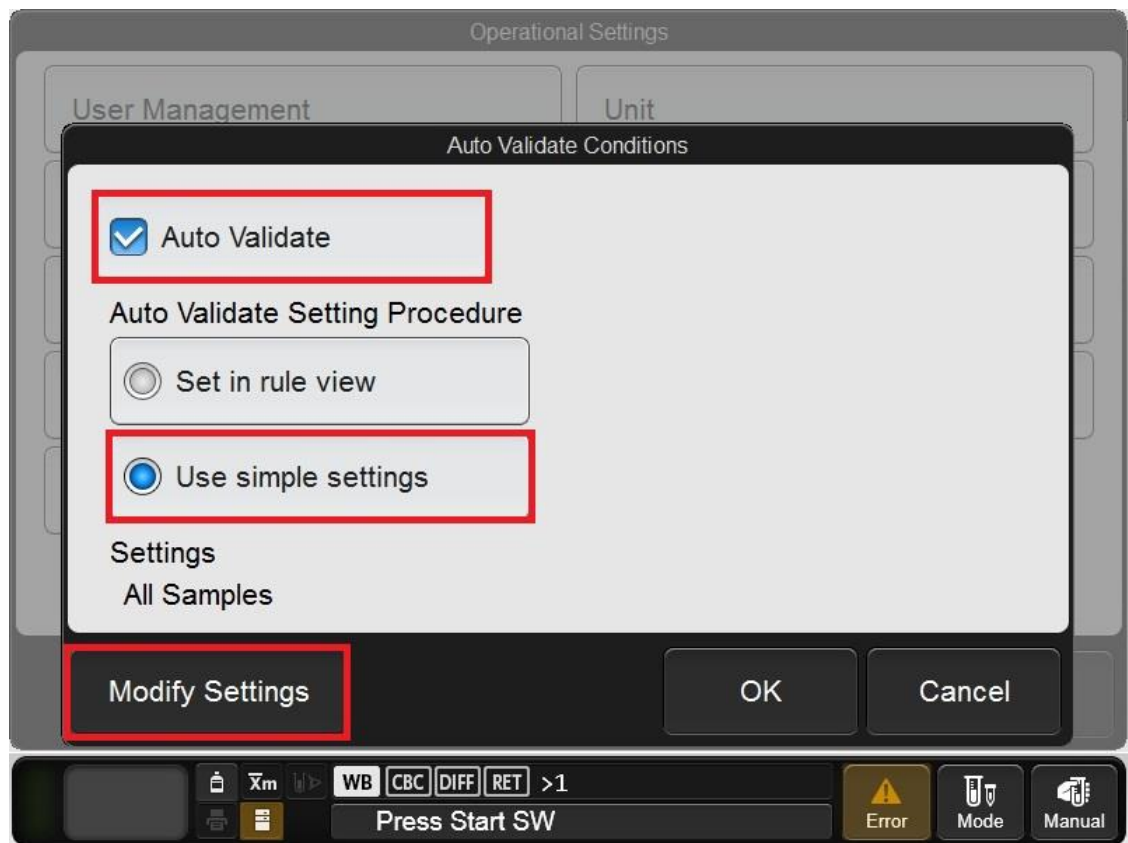


6. Set the settings as shown below
 - a. Select – All Samples
 - b. Click on OK button to apply the setting.

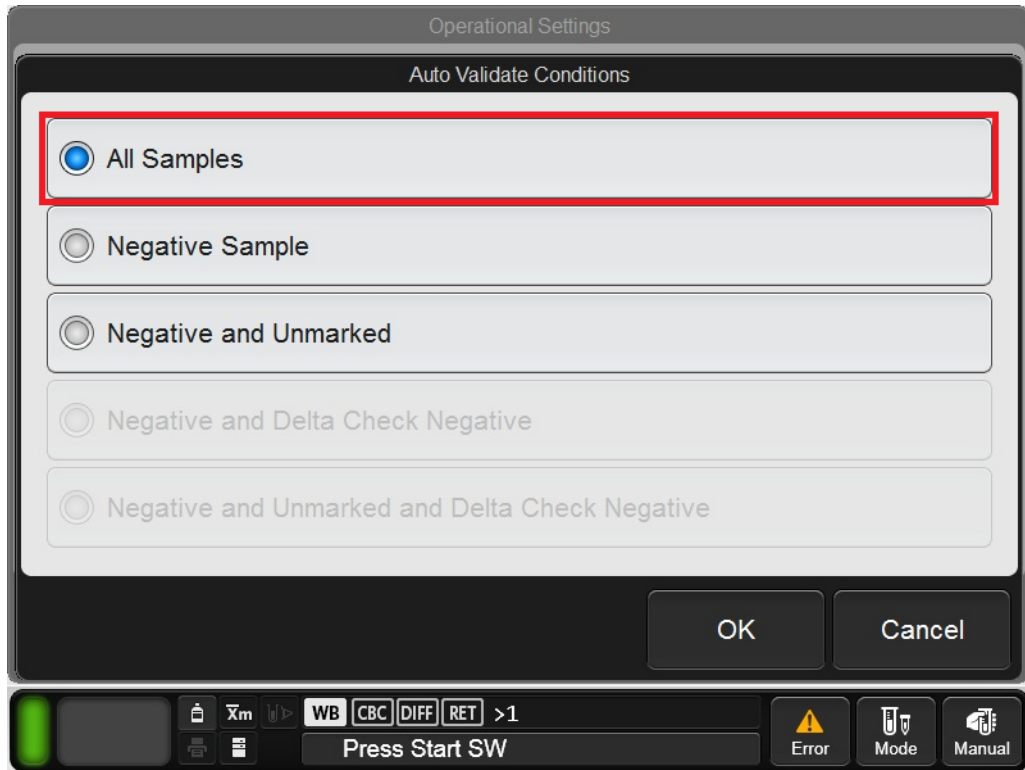


- c. Click on OK button to apply the setting.

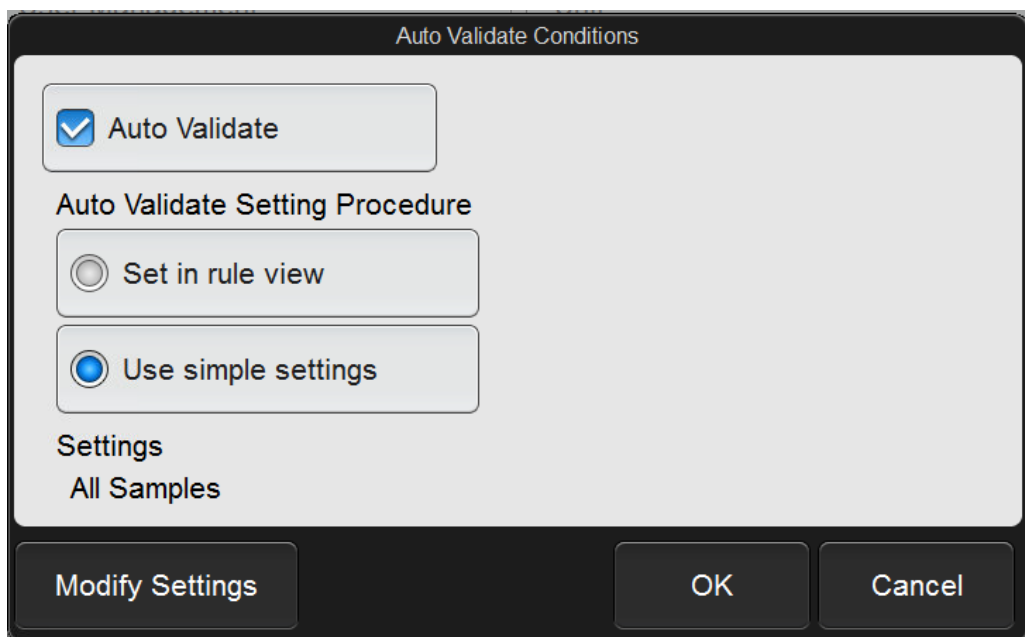


1. Go back to Operational Settings menu,
2. Click on Analysis Ordering

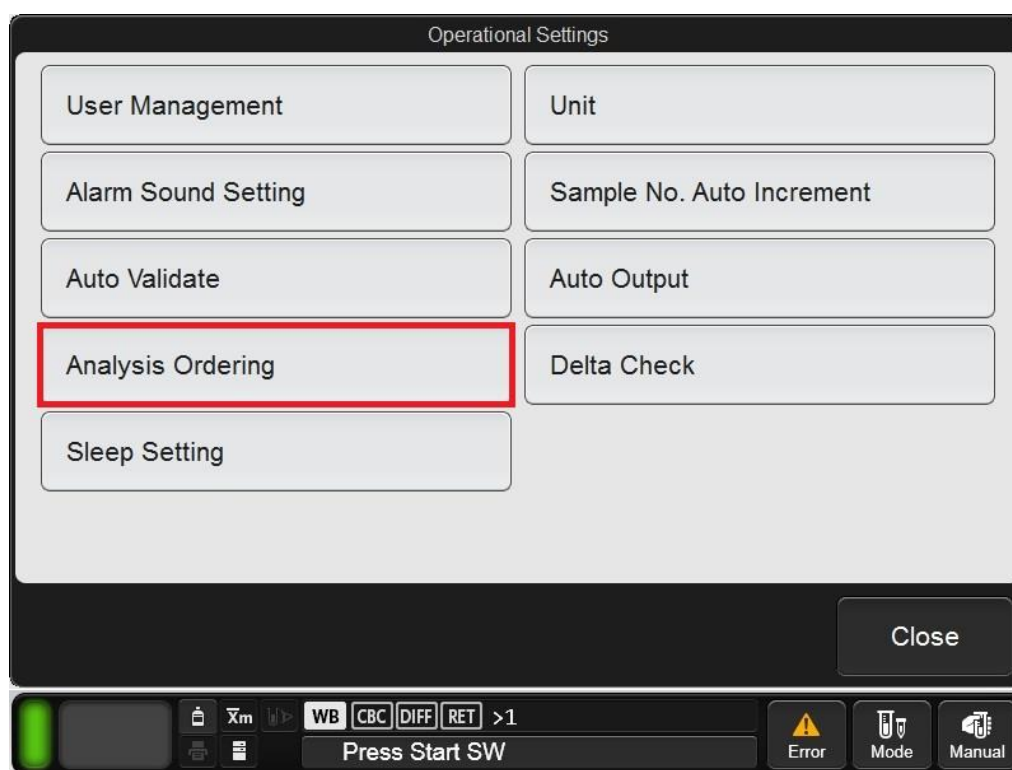


1. Tick – Manual Analysis (If the order key is the sample number).
Note : For XN-550 under Realtime Host Order Setting, Tick – Sample Analysis option.
2. Click on OK button to apply the setting.

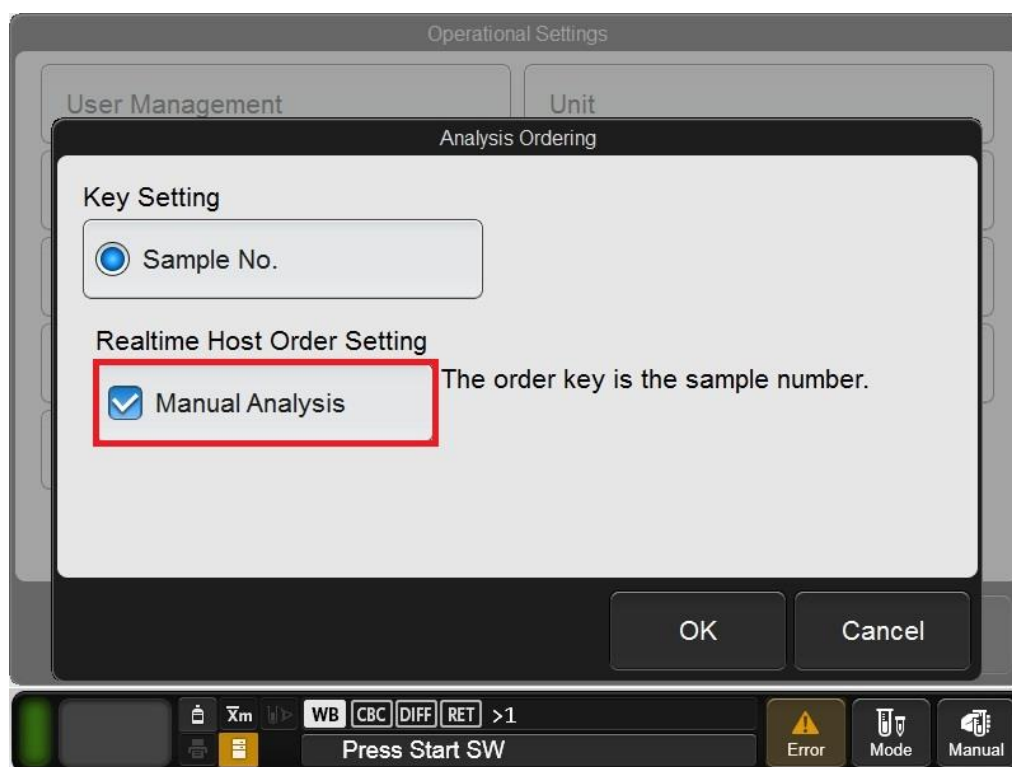

1. Go back to Operational Settings menu,
2. Click on Auto Output

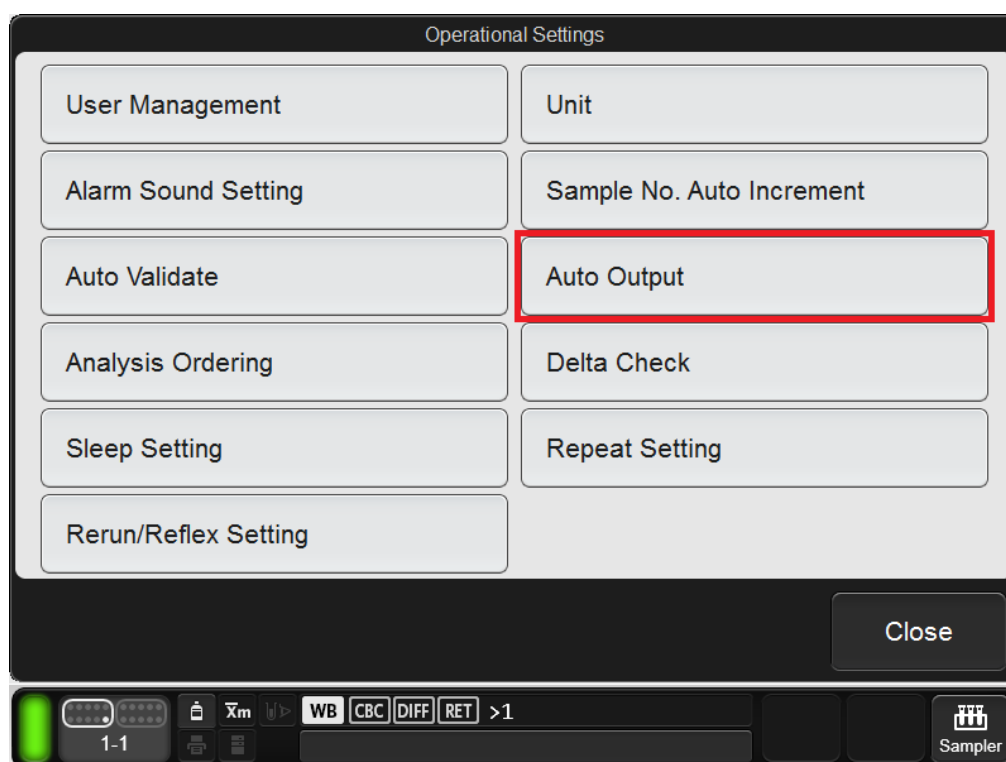


3. Tick – Auto Output
4. Select – Use simple settings

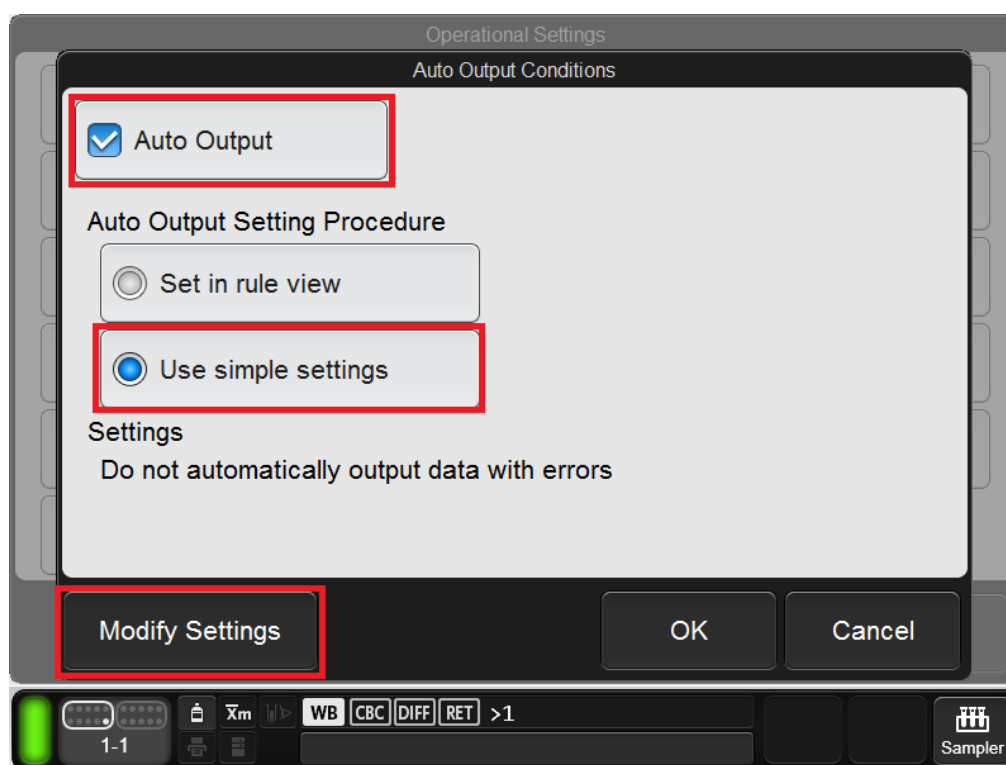

5. Click - Modify Settings button


In the new window below

6. Tick – Do not automatically output data with errors.
7. Auto Output Destination and Output Conditions,
 - a. Tick – HC
 - b. Tick Negative Data – Output
 - c. Tick Diff. Posi. – Output
 - d. Tick Morph. Posi. – Output
 - e. Tick Count Posi. - Output

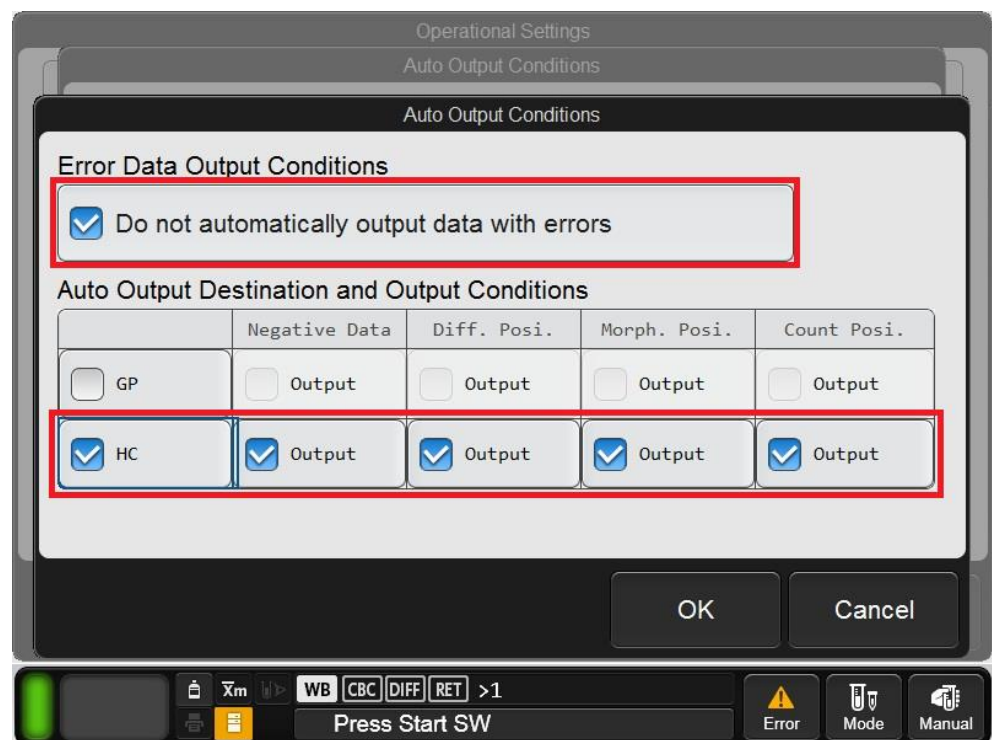


8. Click on OK button to apply the setting.

1. Log in the service account in the IPU software.

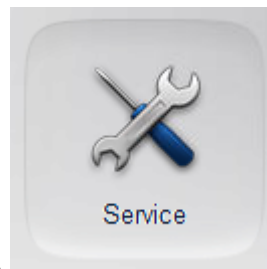


2. Click on Service button.

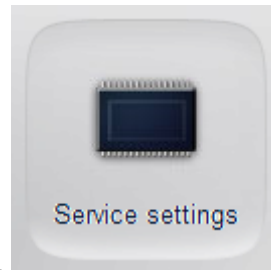

3. Click on Service settings button.

4. Click on HOST

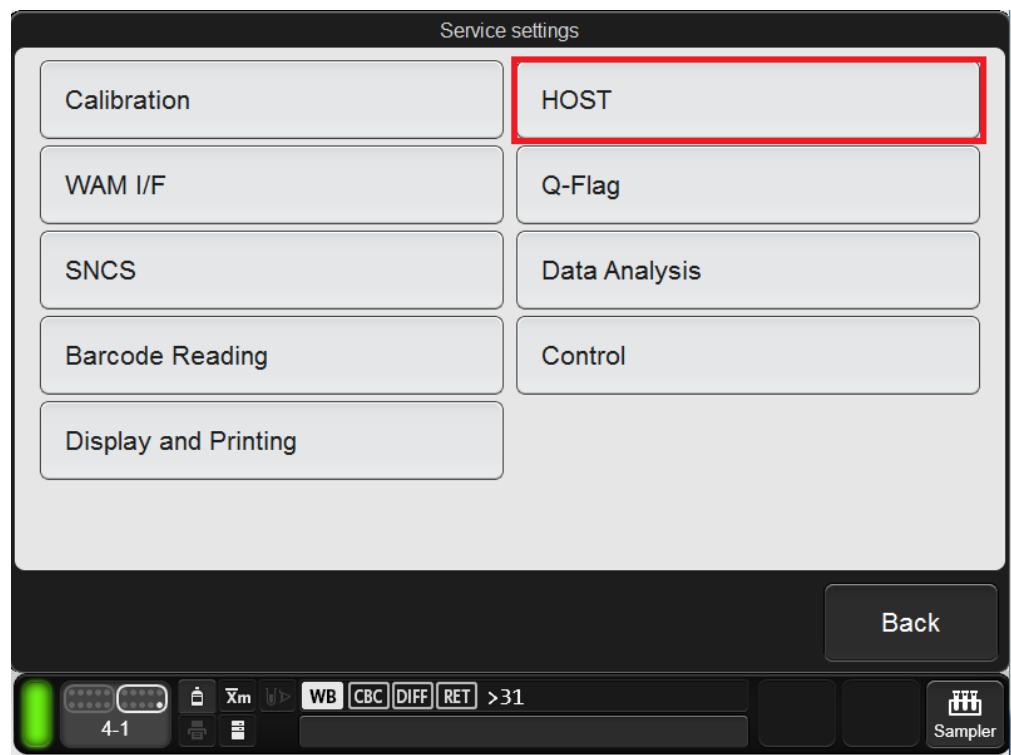

5. Check order key matching to 1. As shown below.

The image shows a 'HOST' configuration window. At the top is a table with two columns: 'Item' and 'Value'. The first row is highlighted in blue. Below the table, there is a section for the selected item with fields for 'Item', 'Description', and 'Value'. The 'Value' field contains the number '1'. At the bottom right are 'OK' and 'Cancel' buttons.

Item	Value
Check order key matching	1
Suppress zeros in sample numbers	0
Delete spaces on right side of sample numbers	0
Sample number padding method when sending an order	0
Sample number padding method for output of analysis results	0
Enable selection of special formats	0

Item: Check order key matching

Description: When communicating analysis information, checks if the value of the order key sent to the host computer matches the value of the order key returned by the host computer. [0: Do Not Check, 1: Check]

Value: 1

Default

OK Cancel

6. Output analysis information to 1. As shown below.

The image shows a 'HOST' configuration window. At the top is a table with two columns: 'Item' and 'Value'. The first row is highlighted in blue. Below the table, there is a section for the selected item with fields for 'Item', 'Description', and 'Value'. The 'Value' field contains the number '1'. At the bottom right are 'OK' and 'Cancel' buttons.

Item	Value
Output analysis information	1
Output research/service parameters	3
Response timeout time	30
Query to host wait time	10
Maximum number of characters in one record (ASTM TCP/IP)	63993
Maximum number of characters in one record (ASTM Serial)	240

Item: Output analysis information

Description: Sets whether analysis information parameters are output when analysis results are output. [0: Do Not Output, 1: Output]

Value: 1

Default

OK Cancel

7. Select “Output research / service parameters” and set value to 3.
This will enable research parameters and service parameters output to host.

The image shows a 'HOST' configuration window. At the top is a table with two columns: 'Item' and 'Value'. The table contains the following items: 'Output analysis information' (1), 'Output research/service parameters' (3), 'Response timeout time' (30), 'Query to host wait time' (10), 'Maximum number of characters in one record (ASTM TCP/IP)' (63993), and 'Maximum number of characters in one record (ASTM Serial)' (240). The 'Output research/service parameters' row is highlighted in blue. To the right of the table are up and down arrow buttons. Below the table, there is a section for the selected item. It has a label 'Item' with the text 'Output research/service parameters'. Below that is a 'Description' box containing the text: 'Sets whether research parameters and service parameters are output when analysis results are output. [0: Do Not Output, 1: Output Research Parameters, 2: Output Research Parameters and Service on IPU Parameters, 3: Output Research Parameters, Service on IPU Parameters, and Service only Host Parameters]'. Below the description is a 'Value' field with the number '3' and a 'Default' button. At the bottom right are 'OK' and 'Cancel' buttons.

Item	Value
Output analysis information	1
Output research/service parameters	3
Response timeout time	30
Query to host wait time	10
Maximum number of characters in one record (ASTM TCP/IP)	63993
Maximum number of characters in one record (ASTM Serial)	240

Item: Output research/service parameters

Description: Sets whether research parameters and service parameters are output when analysis results are output. [0: Do Not Output, 1: Output Research Parameters, 2: Output Research Parameters and Service on IPU Parameters, 3: Output Research Parameters, Service on IPU Parameters, and Service only Host Parameters]

Value: 3

Default

OK Cancel

8. Select "Query to host wait time" and set value to 10.
This will enable the retest judgment to "Query to HOST" for 10 seconds and this will set the wait time from output of the analysis result until the next ordered order query.

"Query host computer for Rerun/Reflex analysis" value should be set to 1(for the above to work)

The image shows a software configuration window titled "HOST". It contains a table with settings and a detailed view for the selected item.

Item	Value
Output analysis information	1
Output research/service parameters	3
Response timeout time	30
Query to host wait time	10
Maximum number of characters in one record (ASTM TCP/IP)	63993
Maximum number of characters in one record (ASTM Serial)	240

Below the table, the "Query to host wait time" item is selected. The "Item" field shows "Query to host wait time". The "Description" field contains the text: "When the retest judgment is 'Query to host', this sets the wait time from output of the analysis result until the next order query (units: sec)". The "Value" field is a text box containing "10". To the right of the value field is a "Default" button. At the bottom right of the dialog are "OK" and "Cancel" buttons.

9. For output of Scattergram and Histogram, set the respective Scattergram/Histogram value to 1. As shown below.

The image shows two screenshots of the 'HOST' configuration window. The top screenshot shows the 'Output WDF(SSC-FSC) Scattergram' item selected in the list, with its value set to 1. The bottom screenshot shows the 'Query host computer for Rerun/Reflex analysis' item selected in the list, with its value set to 1.

HOST

Item	Value
Output WDF(SSC-FSC) Scattergram	1
Output WDF(FSC-SFL) Scattergram	1
Output WDF(FSCW-FSC) Scattergram	1
Output WDF-CBC(FSCW-FSC) Scattergram	1
Output RET(SFL-SSC) Scattergram	1
Output RET(SSC-FSC) Scattergram	1

Item: Output WDF(SSC-FSC) Scattergram

Description: Sets whether WDF(SSC-FSC) Scattergram is output when analysis results are output. (only valid when ASTM protocol is selected)
[0: Do Not Output, 1: Output]

Value: 1

Default

OK **Cancel**

HOST

Item	Value
Output RET(SSC-FSC) Scattergram	1
Output RET(FSCW-FSC) Scattergram	1
Output RBC Histogram with normal range	1
Output PLT Histogram with normal range	1
Query host computer for Rerun/Reflex analysis	1
The method of image data output in the ASTM format	0

Item: Query host computer for Rerun/Reflex analysis

Description: Sets whether the host computer is queried for Rerun/Reflex analysis. [0: Do Not Query, 1:Query, 2: Query Using QueryToHost function in Rule Screen]

Value: 1

Default

OK **Cancel**

10. Select “Query host computer for Rerun/Reflex” value to 1.
This will enable host computer queried for Rerun/Reflex analysis.

The screenshot shows the 'HOST' configuration window. At the top is a table with two columns: 'Item' and 'Value'. The table contains the following items and values:

Item	Value
Output RET(SSC-FSC) Scattergram	1
Output RET(FSCW-FSC) Scattergram	1
Output RBC Histogram with normal range	1
Output PLT Histogram with normal range	1
Query host computer for Rerun/Reflex analysis	1
The method of image data output in the ASTM format	0

Below the table, the 'Item' field is set to 'Query host computer for Rerun/Reflex analysis'. The 'Description' field contains the text: 'Sets whether the host computer is queried for Rerun/Reflex analysis. [0: Do Not Query, 1:Query, 2: Query Using QueryToHost function in Rule Screen]'. The 'Value' field is set to '1'. There is a 'Default' button next to the value field. At the bottom right are 'OK' and 'Cancel' buttons.

11. The method of image data output in the ASTM format” value to 0
(Default folder for PNG path is located at d:\shared)

(Share out this folder on the IPU. Check that the HCLAB WAM computer can access to this shared folder on the IPU.)

The screenshot shows the 'HOST' configuration window. The table at the top is the same as in the previous screenshot, but the 'The method of image data output in the ASTM format' item is now highlighted in blue, indicating it is the selected item. Its value is 0.

Below the table, the 'Item' field is set to 'The method of image data output in the ASTM format'. The 'Description' field contains the text: 'Sets the method of image data output in the ASTM format. (only valid when ASTM protocol is selected) [0: PNG File(Folder Sharing), 1: Raw Data]'. The 'Value' field is set to '0'. There is a 'Default' button next to the value field. At the bottom right are 'OK' and 'Cancel' buttons.

[EOD]