

Document No.	YHLO/ CS -iF3K-20004	Document Version	A2
Application Model	iFlash 1800, iFlash 3000		
Manual Name			
iFlash LIS configuration guide			
Reference Document			
Document No.	Document Name	Version	
Revision Record			
Version	Change Content	Modifier and Date	
A1	First editing	Archy 2018.09.18	
A2	Updating	Allen 2021.10.15	

iFlash LIS configuration guide

A2

Contents

1.1 Foreword.....	4
Chapter I Overview.....	4
1.2 Purpose.....	4
1.3 Scope of document.....	4
Chapter II LIS connection configure.....	4
2.1 The “Communicate Through” is “TCP/IP”	4
2.2 The “Communicate Through” is “Serial Port”	5
Chapter III LIS setting.....	6
3.1 Interface options description.....	6
3.2 Setting option.....	7
3.3 Set Assay Channel No.....	8
3.4 LIS connection.....	9
3.5 Cautious.....	9



Chapter I Overview

1.1 Foreword

With LIS configuration, the instrument can download sample order information from LIS server, or automatically send sample results and quality control results to LIS server in real time or manually in single or batch. LIS server can actively send sample application information to the instrument.

1.2 Purpose

LIS configuration described in this document is used for data exchange between iFlash and laboratory information system (LIS).


1.3 Scope of document

This document is suitable for iFlash 1800 and iFlash 3000

Chapter II LIS connection configure

2.1 The “Communicate Through” is “TCP/IP”

1. Click “Utility --- System Settings --- LIS”

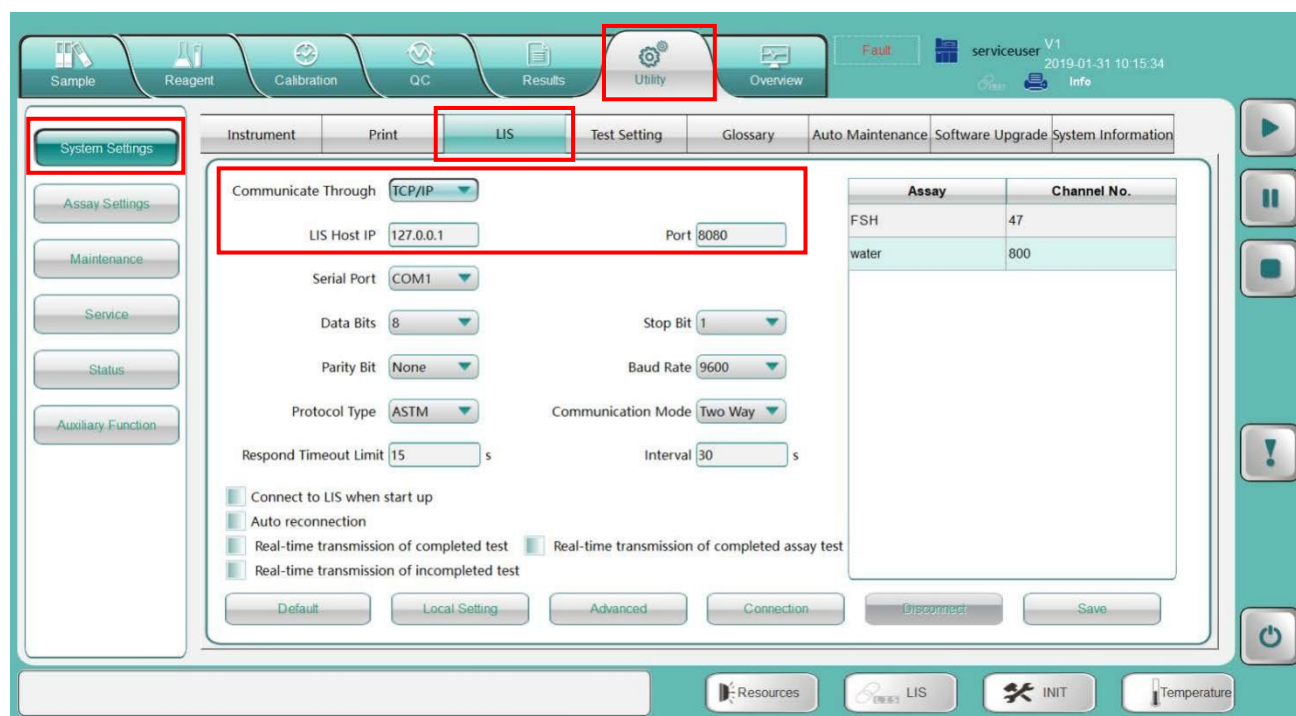


2. Select the “Communicate Through” as "TCP/IP", input the correct “LIS Host IP” and “Port”, click the "save", then click “Connection” to activate the setting.

The IP address of the server is the IP address of the LIS server, and the port is the port number monitored by the LIS server.

2.2 The “Communicate Through” is “Serial Port”


1. Click “Utility --- System Settings --- LIS”.

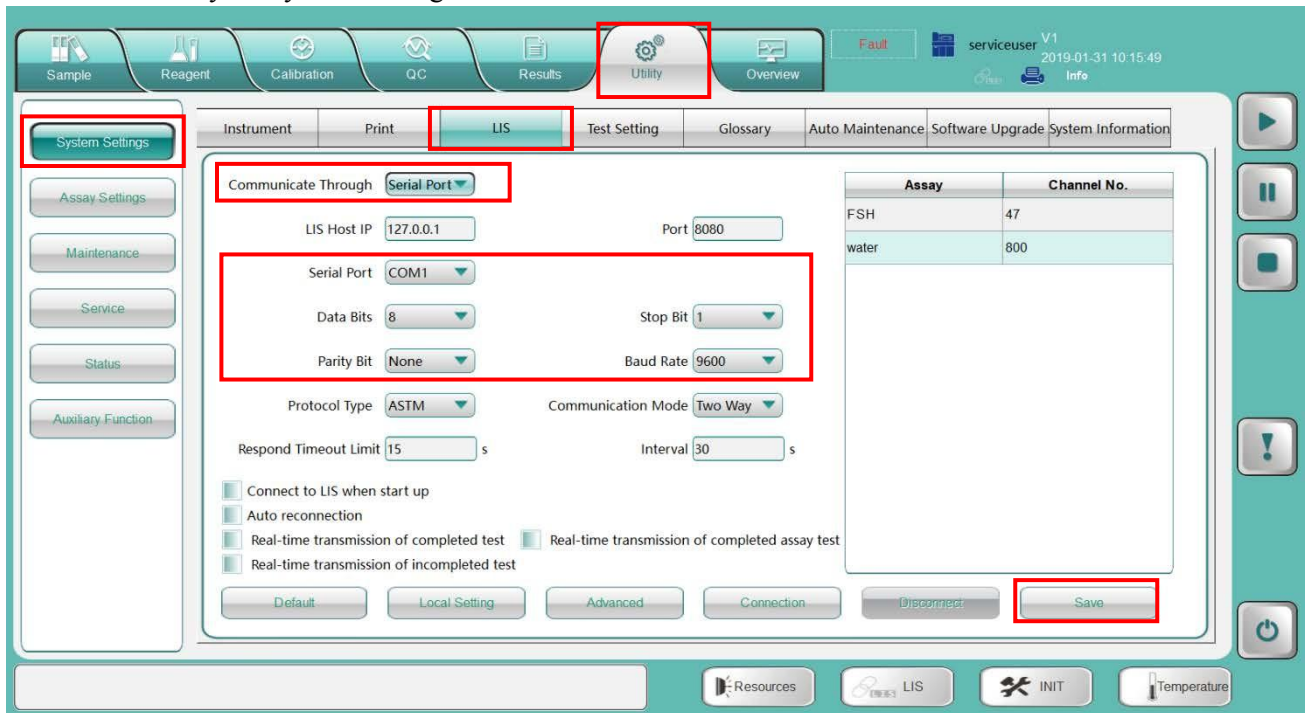

2. Select the “Communicate Through” as "Serial Port", select the correct “Serial port”, “Data Bits”, “Stop Bit”, “Parity Bit” and “Baud Rate” click the "save", then click “Connection” to activate the setting.

Chapter III LIS setting

3.1 Interface options description

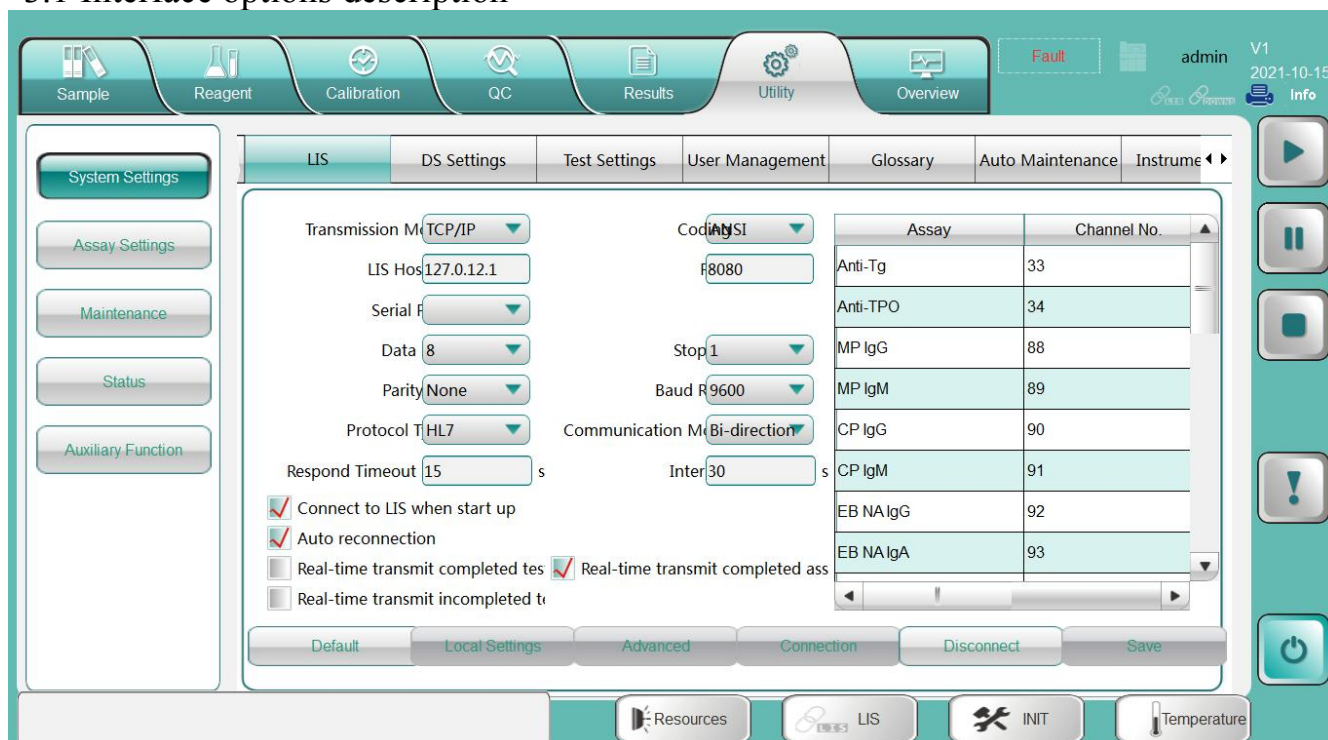

1. [Code] Currently, two encoding types [ANSI] and [UTF-8] are supported, the default is [ANSI], and the encoding method should be consulted by LIS engineers;
2. [Protocol] : select the corresponding protocol according to LIS protocol of LIS engineers. Currently, the software supports two protocol modes: [HL7] and [ASTM], and the main protocol mode of domestic customers is [HL7].
3. [Communication Model] : Select [Bidirectional] or [one-way] according to the actual type.
 [Bidirectional] : LIS server can actively send the application information of samples to the instrument, and the test results can be automatically transmitted to LIS server after the instrument test is completed.
 [One-way] : The device cannot download test information from LIS server, so it supports automatic or manual transmission of test results to LIS server after manually applying samples or completing quality control tests;
4. [Respond timeout] : Set the connection timeout limit between the instrument and LIS, the range is [10-60s], and the default setting is 15s;
5. [Interval] : Enter the interval for automatic reconnection of the system after LIS connection is disconnected, which ranges from [30 to 60s] and is set to 30s by default.
6. [Connect to LIS when start up] : After this item is selected, the system will actively establish connection with LIS server upon startup;
7. [Auto reconnection] : After this item is selected, the system will check LIS connection status at the set [interval] . If LIS is disconnected, and LIS is automatically reconnected; if the connection is normal, LIS is not reconnected;
8. [Real-time transmit completed test] : After this item is selected, when the state of a sample changes from [testing] to [completed], the system will automatically send the test results to LIS. This function is only applicable to the sample test results of the day, and the historical sample results cannot be automatically sent to LIS.
9. [Real-time transmit incompleted test] : After this item is selected, when the state of a sample changes

from

[testing] to [unfinished], the system will automatically send the test results to LIS. This function is only applicable to the sample test results of the day, and the historical sample results cannot be automatically sent to LIS.

10. [Real-time transmit completed assay] : After this item is selected, when the status of the test project changes from [testing] to [completed], the system will automatically send the test results to LIS. This function is only applicable to the sample test results of the day, and the historical sample results cannot be automatically sent to LIS.

11. [Automatic numbering of unnumbered samples (starting from 100001)] : After this item is selected, when LIS downloads the sample information without sample number, the samples will be numbered from 100001 automatically;

12. [LIS download automatic retest] : Only visible under ServiceUser and YHLO permissions (not visible under admin permissions). After selecting this item, the test results of the sample will be deleted at LIS terminal and retested. The instrument will be retested after receiving the test application issued by LIS; You are advised not to select this option.

3.2 Setting option.

1) Data transmission according to sample, and the setting method is shown as follows:

LIS Settings	
municate Through	TCP/IP
LIS Host IP	127.0.12.1
Serial Port	
Data Bits	8
Parity Bit	None
Protocol Type	HL7
Send Timeout Limit	15 s
Code	ANSI
Port	8080
LIS Instrument No.	iFlash3000
Stop Bit	1
Baud Rate	9600
Communication Mode	Two-Way
Interval	30 s
☒ Connect to LIS when start up ☒ Auto Reconnection ☒ Real-time upload completed test ☒ Real-time upload incompleted te	
☐ Unnumbered samples auto num ☐ Real-time upload assay results	

Buttons: Default, Local Setting, Connection, Disconnect, Save

2) Data transmission according to assay, the setting method is as follows:


3.3 Set Assay Channel No.

Set the corresponding relationship between the assay and LIS server, and LIS server displays the corresponding channel number.

- 1) Add Assay: **The channel number of the assay is same with the assay ID number.** The channel number of the assay can be displayed by adding the item.
- 2) Add methodology : [Utility] --[Assay Settings] -- [Assay] -- [Add]

3.4 LIS connection

1. Click “Connection”, LIS connection status is successful, as shown below:

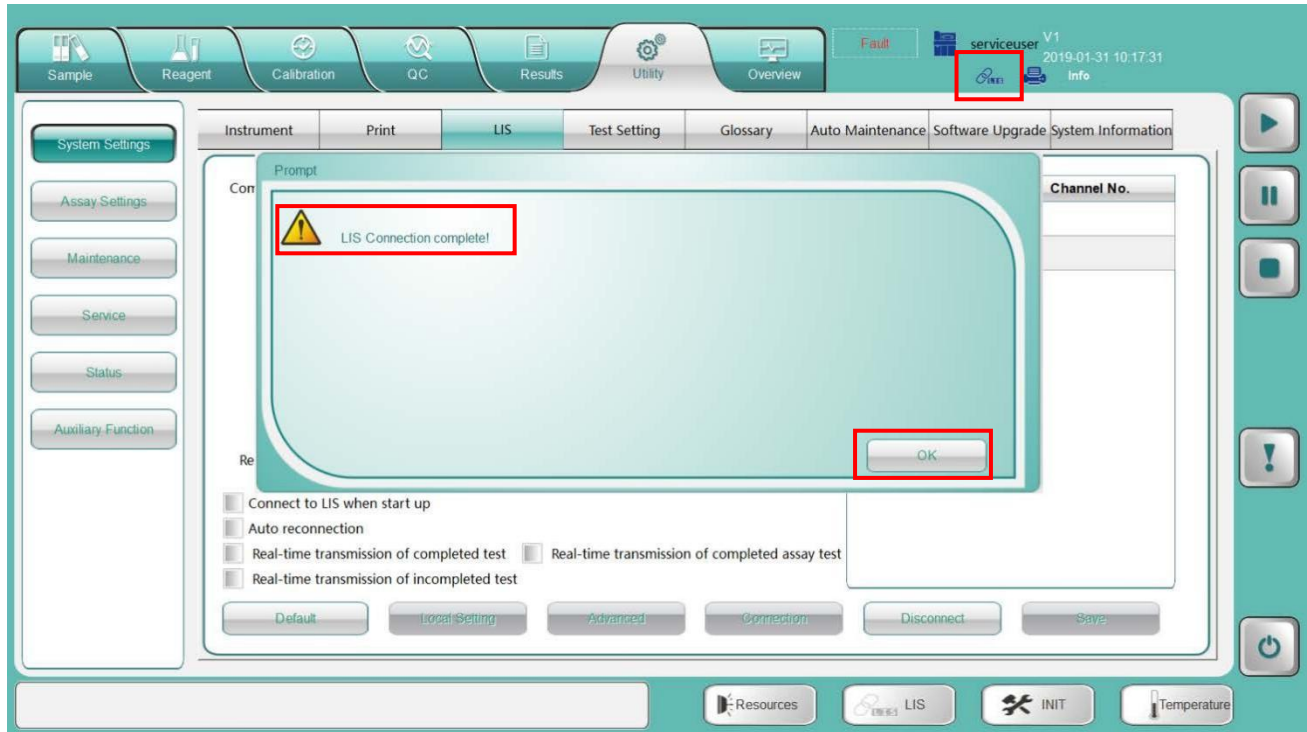

2. Click “OK” button to close the dialog interface

3.5 Cautious

1. Disconnect LIS connection first before modify LIS setting
2. Click the “SAVE” button after modify finish, then click “Connection” to activate the setting.