



TECHNICAL BULLETIN

Ref: TB 062

Author: M. Legrand

Validation: J. Stawicki

For:

IMMEDIATE ACTION	☒
ACTION	☐
UPDATE	☐
INFORMATION	☐


	Section
Model: MYHTIC 22 AL	☛ 1
Information: Bug causing impossibility to restart the instrument.	☛ 2
Action: Do not set "DELAY ms" parameter lower than 100 millisecond.	☛ 3
Send by ORPHE SA to: All distributors	
Contacts: Contact ORPHEE SA for further information.	
Information feedback asked by ORPHEE SA To inform ORPHEE SA of the incidents relating to this document.	

1 MODEL

All Mythic 22 AL instruments.

2 INFORMATION

Orpheediscovered a bug which causes an impossibility to restart the instrument. Actually, when the “**DELAY ms**” parameter from “**COMMUNICATION**” screen is set with a lower value than 100, it is impossible to restart the instrument after shutdown.

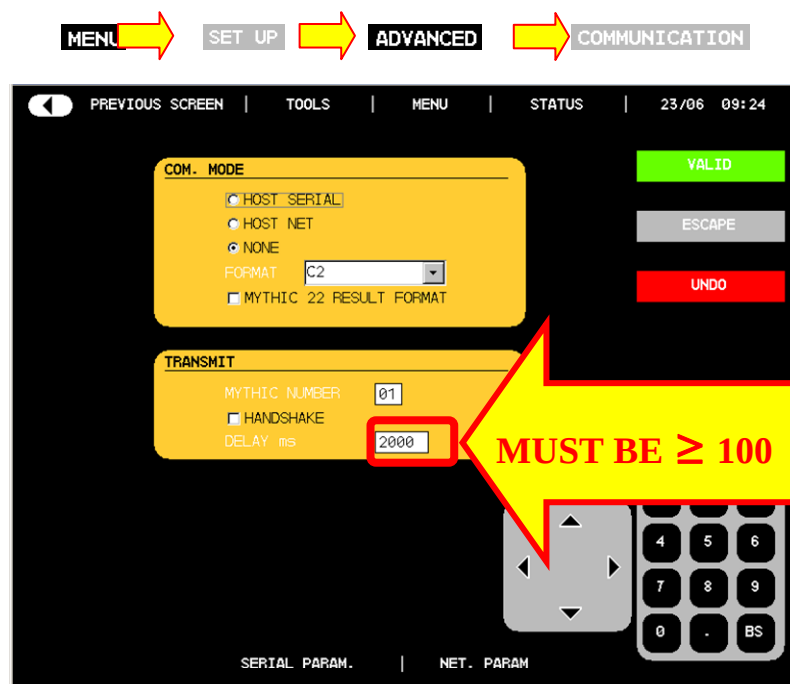

3 IMMEDIATE ACTION

In order to avoid this unpleasant scenario, it is strongly recommended to inform all distributors to always set this parameter with a value **higher or equal to 100**. The bug will be corrected with a future software release.