

Enriching the quality of life for people everywhere.

The easy way to test the iQ computer's LIS port

Greg Loa

October 17, 2007



Iris Diagnostics

A Division of IRIS International, Inc.

The easy way to test the iQ computer's LIS port

- What is the purpose?
- What is needed?
- How to do it?

What is the purpose?

- The LIS port is an integral part of the iQ computer and must be functioning.
- The LIS port could be on, but the iQ application may not transmit results.
- The iQ results can't be transmitted via network.

What is needed?

1. Your Laptop.
2. One DB9 FF Null Modem cable.
(P/N: 250-3109)
3. Few minutes of your time.



How to do it?

1. Connect one end of the cable to the LIS port of the iQ Computer.



2. Connect the other end of the cable to the serial port of your laptop or computer.



3. Configure the Hyper Terminal in your laptop or computer by doing the following:

Left click “Start”


Choose:

Programs


Accessories

Communications

Left click on “Hyper Terminal”

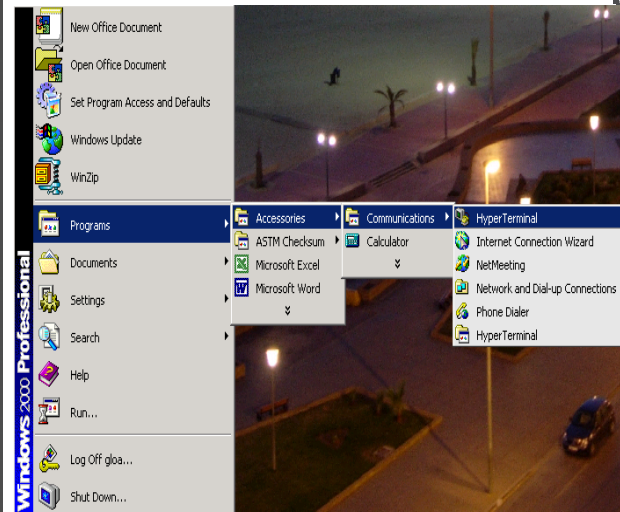


If you have Window XP, click “No” on “Default Telnet Program?” and you will see a screen shown on this page.

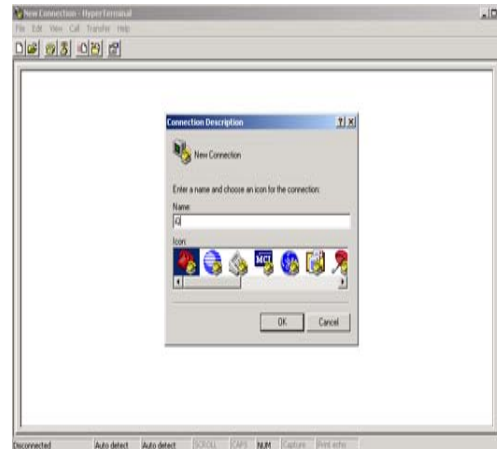


Name the new connection, such as iQ

Click “OK”.

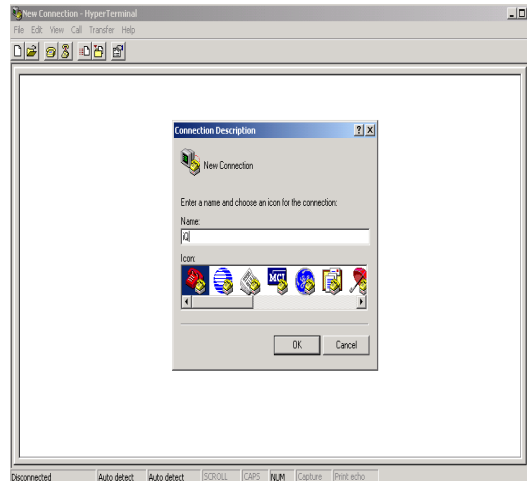


Choose connect using: COM1 and
click “OK”.

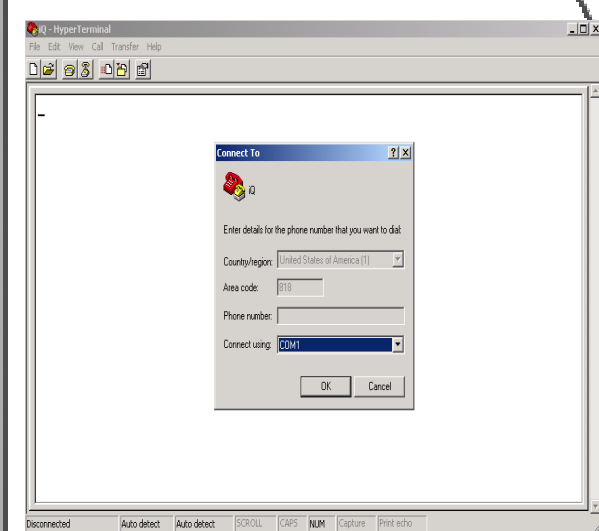

Configure COM1:


Bits per second 9600
Data bits 8
Parity None
Stop bits 1
Flow control None

Then click “OK”.

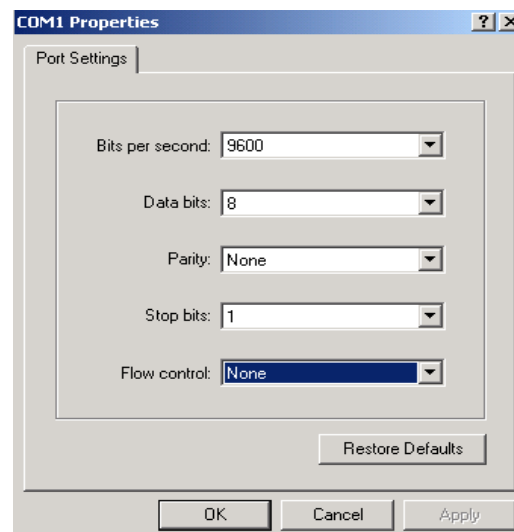

The Hyper Terminal in your laptop is now ready to receive transmission from the iQ.

Go to the iQ screen and follow the instruction in the next slide.

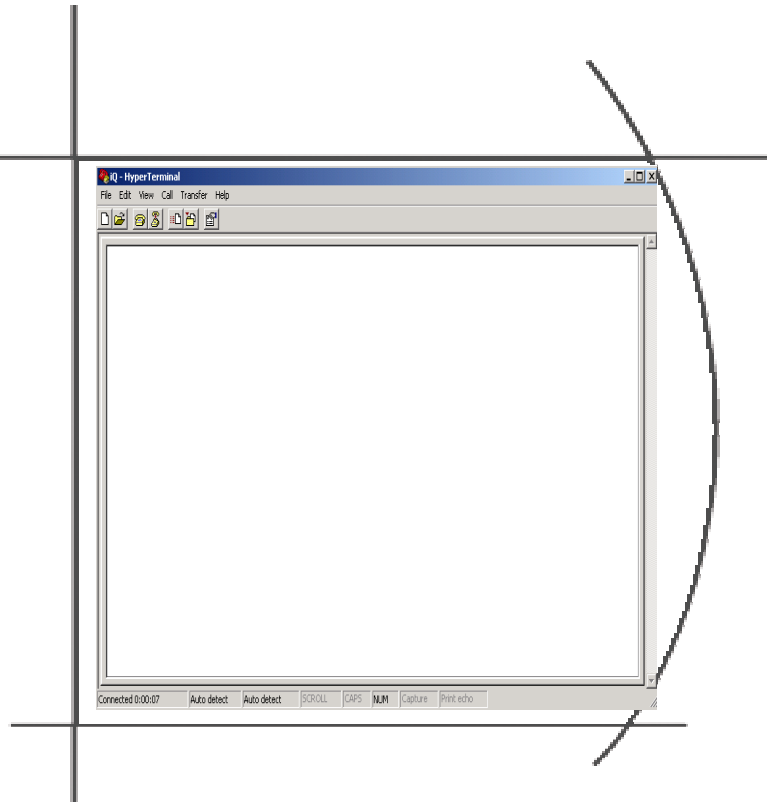


From the iQ (v5.0) Main menu:

Logon as Tech or Manager (Do not use [IRIS])

Click “**Work List**”

Click “**Search**”

Enable “**Show released specimen only**” and click “OK”. You will see the “Found List” as shown on this page.

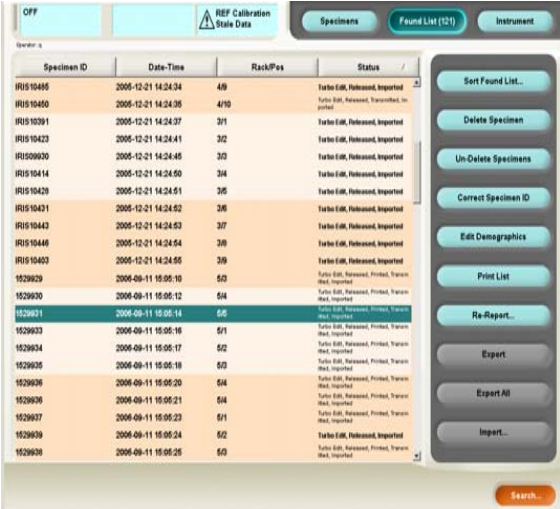

Specimen ID	Date-Time	Rack/Pos	Status
IRIS10485	2005-12-21 14:24:34	4/9	Turbo I-Off, Released, Imported
IRIS10486	2005-12-21 14:24:35	4/10	Turbo I-Off, Released, Transmitted, Imported
IRIS10391	2005-12-21 14:24:37	3/1	Turbo I-Off, Released, Imported
IRIS10423	2005-12-21 14:24:41	3/2	Turbo I-Off, Released, Imported
IRIS09030	2005-12-21 14:24:45	3/3	Turbo I-Off, Released, Imported
IRIS10414	2005-12-21 14:24:50	3/4	Turbo I-Off, Released, Imported
IRIS10428	2005-12-21 14:24:51	3/5	Turbo I-Off, Released, Imported
IRIS10431	2005-12-21 14:24:52	3/6	Turbo I-Off, Released, Imported
IRIS10443	2005-12-21 14:24:53	3/7	Turbo I-Off, Released, Imported
IRIS10446	2005-12-21 14:24:54	3/8	Turbo I-Off, Released, Imported
IRIS10403	2005-12-21 14:24:55	3/9	Turbo I-Off, Released, Imported
1529629	2006-09-11 15:05:10	5/3	Turbo I-Off, Released, Printed, Turbomax, Imported
1529630	2006-09-11 15:05:12	5/4	Turbo I-Off, Released, Printed, Turbomax, Imported
1529631	2006-09-11 15:05:14	5/5	Turbo I-Off, Released, Printed, Turbomax, Imported
1529632	2006-09-11 15:05:16	5/1	Turbo I-Off, Released, Printed, Turbomax, Imported
1529634	2006-09-11 15:05:17	5/2	Turbo I-Off, Released, Printed, Turbomax, Imported
1529635	2006-09-11 15:05:18	5/3	Turbo I-Off, Released, Printed, Turbomax, Imported
1529636	2006-09-11 15:05:20	5/4	Turbo I-Off, Released, Printed, Turbomax, Imported
1529636	2006-09-11 15:05:21	5/4	Turbo I-Off, Released, Printed, Turbomax, Imported
1529637	2006-09-11 15:05:23	5/1	Turbo I-Off, Released, Printed, Turbomax, Imported
1529639	2006-09-11 15:05:24	5/2	Turbo I-Off, Released, Imported
1529638	2006-09-11 15:05:25	5/3	Turbo I-Off, Released, Printed, Turbomax, Imported

From the iQ (v5.0) Found List:

Highlight a result


Click **"Re-Report"**

Disable **"Screen"**

Enable **"LIS"**

Click **"OK"**

Go to your laptop or
computer immediately!!

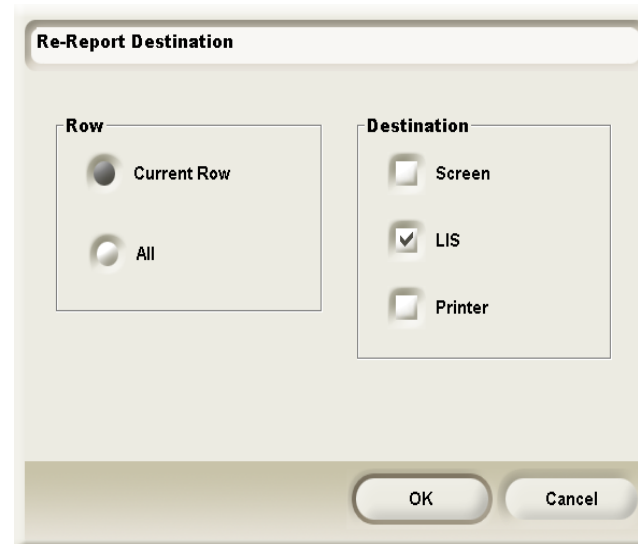

The image shows a 'Re-Report Destination' dialog box. It has two main sections: 'Row' and 'Destination'. The 'Row' section has two radio buttons: 'Current Row' (selected) and 'All'. The 'Destination' section has three checkboxes: 'Screen' (unchecked), 'LIS' (checked), and 'Printer' (unchecked). At the bottom right, there are 'OK' and 'Cancel' buttons.

Re-Report Destination	
Row	Destination
☒ Current Row	☐ Screen
☐ All	☒ LIS
	☐ Printer
OK Cancel	


While looking at the Hyper Terminal screen in your laptop, click either **both “Ctrl” and “F” keys at the same time**

or

press and hold “Ctrl” key and press “F” key

one frame at a time.

You will see iQ result transmission as shown in the next slide.

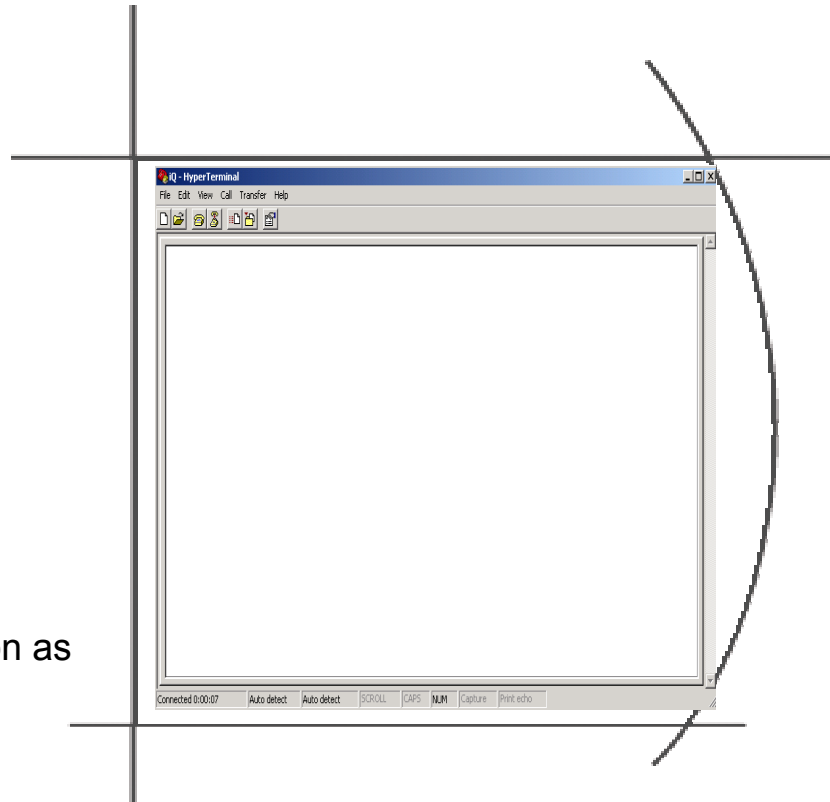


The slide on this page shows how the iQ result transmission looks like.

You are not required to understand, however, it is important for you to see that the iQ computer is transmitting results to the LIS.

```
0061006C002200200041004600300022003000220020004E0052003000220046004500570022003E
001E2
005B006E006F006E0065005D003C002F00410052003E003C004100520020004B0065007900300022
0043004100430042002200200053004FE003000220043004100430042002200004C004FE00300022
00430061006C006300690075006D000200043006100720062006F006E0061007400650062000430072
001E7
0179007300740061006C002200200041004600300022003000220020004E00520030002200460045
00570022003E005B006E006F006E0065005D003C002F00410052003E003C004100520020004B0065
0079003000220055005200490043002200200053004FE00300022005500520049004300220020004C
001E4
024E003000220055007200650063002000410063006900640020004300720079007300740061006C
002200200041004600300022003000220020004E0052003000220046004500570022003E005B006E
006F006E0065005D003C002F00410052003E003C004100520020004B0065007900300022004C0045
001E6
03550043002200200053004FE00300022004C00450055004300220020004C004FE00300022004C0065
00750063006C006E006500220004300720079007300740061006C0022002000410046003000220030
00220020004E0052003000220055004FE00530022003E0065006E006F006E0065005D003C002F0041
001E9
0452003E003C004100520020004B006500790030002200430059005300654002200200053004FE0030
00220043005900530065400220020004C004FE0030002200430079007300740061006C0065006200043
00720079007300740061006C002200200041004600300022003000220020004E00520030002200550
001E0
054F00530022003E005B006E006F006E0065005D003C002F00410052003E003C004100520020004B
0065007900300022005400590052004FE002200200053004FE00300022005400590052004FE00220020
004C004FE00300022005400790072006F00730069006E00650020004300720079007300740061006C
001E0
062200200041004600300022003000220020004E00520030002200550004FE00530022003E005B006E
006F006E0065005D003C002F00410052003E003C004100520020004B00650079003000220041004D
004FE0052002200200053004FE003000220041004D004FE005200220020004C004FE003000220041006D
001E1
076F00720070068006F007500730020004300720079007300740061006C00220020004100460030
0022003000220020004E0052003000220046004500570022003E005B006E006F006E0065005D003C
002F00410052003E003C004100520020004B00650079003000220055004FE0043004C002200200053
001E3
004E003000220055004FE0043004C00220020004C004FE003000220055006FE0063006C006100730073
00690066006C00650064002200200041004600300022003000220020004E00520030002200390039
0059003900590059005900590020002F0075004C0022003E005B006E006F006E0065005D003C002F
001E6
01410052003E003C002F00410043003E003C004100520056003E00340032003C002F004100520056
003E003C00410055004C002F003E003C00530079007300430066006E0020004D006C00630072006F
00500072006500730065006E0074003000220054002200200041005800500072006500730065006E
001E51
0274003000220054002200200069004300680065006D00320030003000500072006500730065006E
00740030002200460022002F003E003C002F00530041003E000D000000#63
```


When done, close your Hyper Terminal session by clicking the “x” in the right upper corner of the Hyper Terminal screen and click “Yes” to disconnect now.

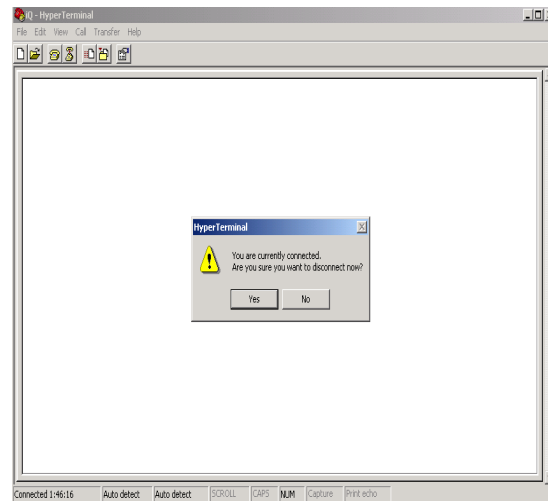

If this connection is the first connection you created, then you will see the screen shown on this page.

Click “Yes” to save the session and Disconnect your cable.

